

HEADRESTS

Tecnodent manual Headrests are available in 4 different models:

Elle Fisso, Elle/Elle2 Mobile, UNI and Light.

Every headrest was made for meeting every operator's need.

Elle/Elle2

ELLE/Elle2 Fisso has only the movement along its vertical axis. ELLE/Elle2 Mobile is an anatomical headrest, characterized by movement of inclination of the upper part (double articulation) for a correct cervical support.



Elle / Elle 2



C 95



C 2002

Magnetic pillows

Elle/Elle2 type, both fix and mobile version, can be combined with magnetic pillow C95 or C2002. Either pillow can be applied directly on the backrest, for ensure a better comfort, independently of the patient's height.



UNI



Light

UNI and Light

The UNI and Light headrests are multi-articulated and allows positioning of the headrest along 3 axis.

Linda3 ORL chair was created to perfectly satisfy professionals in the audiology, otology, ophthalmology, dental surgery, general practice, as well as cosmetic field.

Linda3 ORL is equipped with an electromechanical up&down system and tilting control of the backrest and a backrest-leg rest synchronism.

Linda3 ORL has a rotating base, with rotation angle up to 325°, anatomical and seamless upholsteries with lateral supporting wings, Trendelenburg movement, armrests, memorizable operating positions and anatomical headrest.

All functions are electronically controlled by infra-red device and/or high precision pedal control, so that the doctor is free to operate in total comfort and hygiene with hands free.

Linda3 ORL combines functionality, ergonomics, safety, quality with a smart and reassuring design.

TECHNICAL FEATURES


| | | |
|---|----|-----------|
| Maximum height | mm | 840 |
| Minimum height | mm | 510 |
| Rotation angle | ° | 325 |
| Width | mm | 660 |
| Length | mm | 1900 |
| Maximum lifting weight (patient + unit) | kg | 135 + 75 |
| Gross weight | kg | 145 |
| Net weight | kg | 115 |
| Supply power | V | 230 (115) |
| Frequence | Hz | 50 (60) |
| Insulation class | | I |
| Parts type | | B |

TECNODENT
INNOVATIVE MEDICAL SOLUTIONS

Tecnodent S.r.l. Via 63° Brigata Bolero, 24 – 40033 Casalecchio di Reno (BO), Italy
Tel. +39 051 6131143 – Fax +39 051 575402 – info@tecnodent.com – www.tecnodent.com

TECNODENT
INNOVATIVE MEDICAL SOLUTIONS



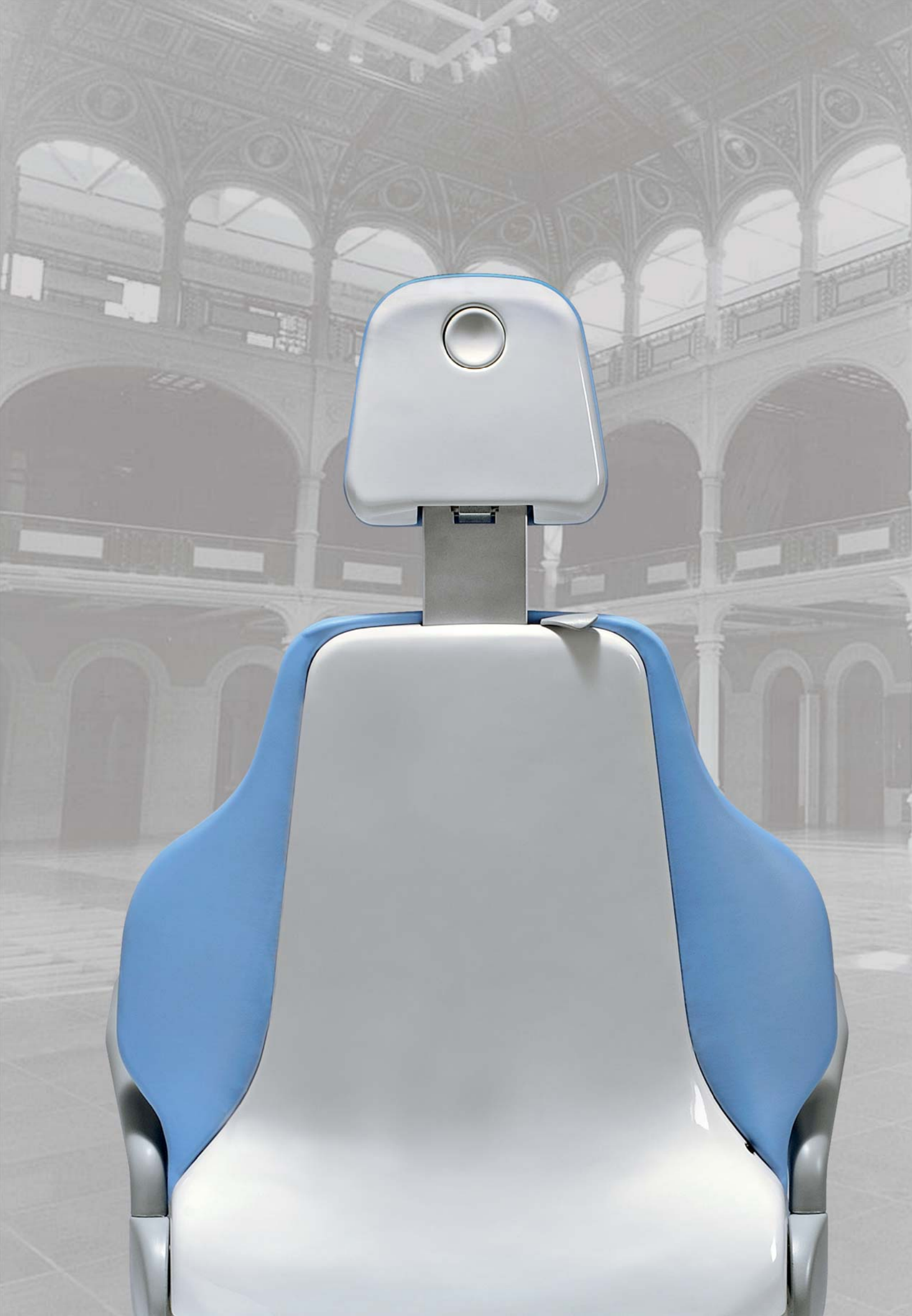

Linda3 ORL

Linda3 ORL guarantees:

- Electromechanical up&down and backrest tilting movements, backrest-leg rest synchronism to reach the Trendelenburg position;
- ELLE fix headrest, with long-section sliding;
- Longitudinal rotating base with rotation angle 325°;
- Anatomical backrest with lateral armrest support wings for the highest patient's comfort;
- Complete symmetry for left- and right-handed professionals;
- Infra-red control;
- Anti-crushing safety stop onto the backrest and leg rest, for an emergency stop of the chair;
- Automatic memorized return to "0" position at the end of the treatment;
- Automatic memory of the last operating position;
- 3 Programmable operating positions;
- Tilting leg rest, inclination up to 90° for get on and off the chair in an easier way;
- Minimum height 50 cm, for the easiest access to the chair, also for elder or disabled persons.



Infrared control

Functioning remote control for chair movements and memorized positions



Mobile pedal control "Easy Touch" (OPTIONAL)

Easy and functioning pedal control, for the highest freedom and hygiene for the doctor during the treatment.

Right and/or left armrest

Rotating armrests help reach comfort and get on and off the chair



Rotating base

The chair can rotate up to 325° along its vertical axis and is provided with rotation lock/unlock.